

# Defined cGMP animal-origin-free medium supports PSC during, and promotes single-cell clonal expansion after, microfluidic sorting




Korina Eribez<sup>2</sup>, Adonary Munoz<sup>2</sup>, Gerardo Aguirre Castillo<sup>1</sup>, Jerome Karpiak<sup>1</sup>
1 Defined Bioscience Inc. 6404 Nancy Ridge Dr, San Diego, CA 92121, 2 Nanocellect Biomedical, 9525 Towne Centre Dr, San Diego, CA 92121

## Introduction

PSC fate is modulated by environmental cues. Dissociation, followed by deposition into a single-PSC microenvironment to establish a cell line, challenges clonal propagation. convention aggregate Subsequent inefficiency of clonal expansion is a barrier to gene editing workflows for basic and translational research. Further, successful translation of stem cell-based therapies requires tracking the cell product and the components in contact with the product throughout the manufacturing process. To facilitate this, the recent launch of our cGMP-compliant, fully defined, animal-origin-free complete HiDef-B8 PSC maintenance and expansion medium kit is in high demand. HiDef-B8-cGMP is an optimized PSC maintenance medium manufactured under relevant cGMPs and was developed with traceable, animal-free raw materials. Here we describe a protocol for robust single-cell cloning using Nanocellect's WOLF G2 microfluidicsbased cell sorter in combination with HiDef-B8cGMP complete medium used as both sheath fluid, as well as PSC maintenance medium.


# **WOLF® G2 Microfluidic Sorter**




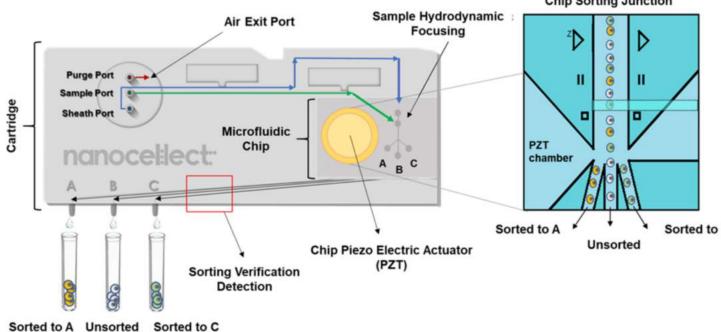

## Results

#### **PSC** maintenance and sorting

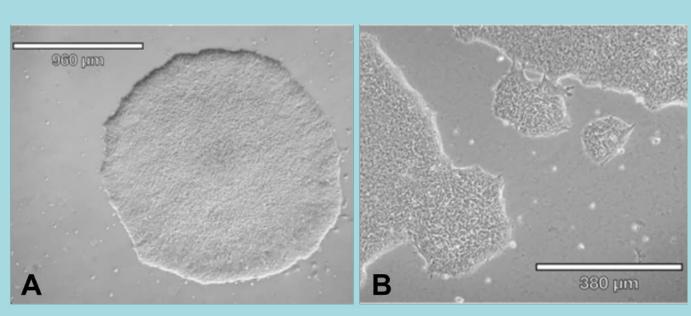

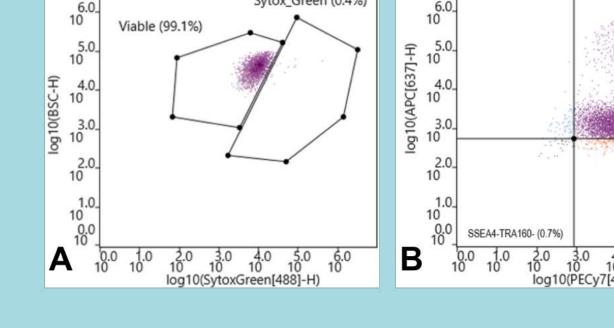

cultured on cGMP-compliant CTS™ Defined VTN-N substrate in Bioscience's DMEM/F-12 formulation TrypLE™ with L-alanyl-L-glutamine and HiDef-**B8-cGMP:** (A) Representative image of a healthy compact colony.

(B) Higher mag representative image well defined edges, and nucleus:cytoplasm ratio.

Figure 2. Morphology of PSCs Figure 3. Single-cell suspensions of PSCs in complete HiDef-B8-cGMP medium with 1% cGMP-compliant Select were gated for viability and dual expression of pluripotency markers TRA-1-60-R and SSEA-4 using WOLFViewer software: Bulk flow cytometry (after WOLF N1 of colonies depicting distinct borders, single-cell deposition) demonstrates large parent sorting populations with (A) 99% viability and (B) 93% dual expression.

#### Initial adherence and clonal outgrowth post-sort

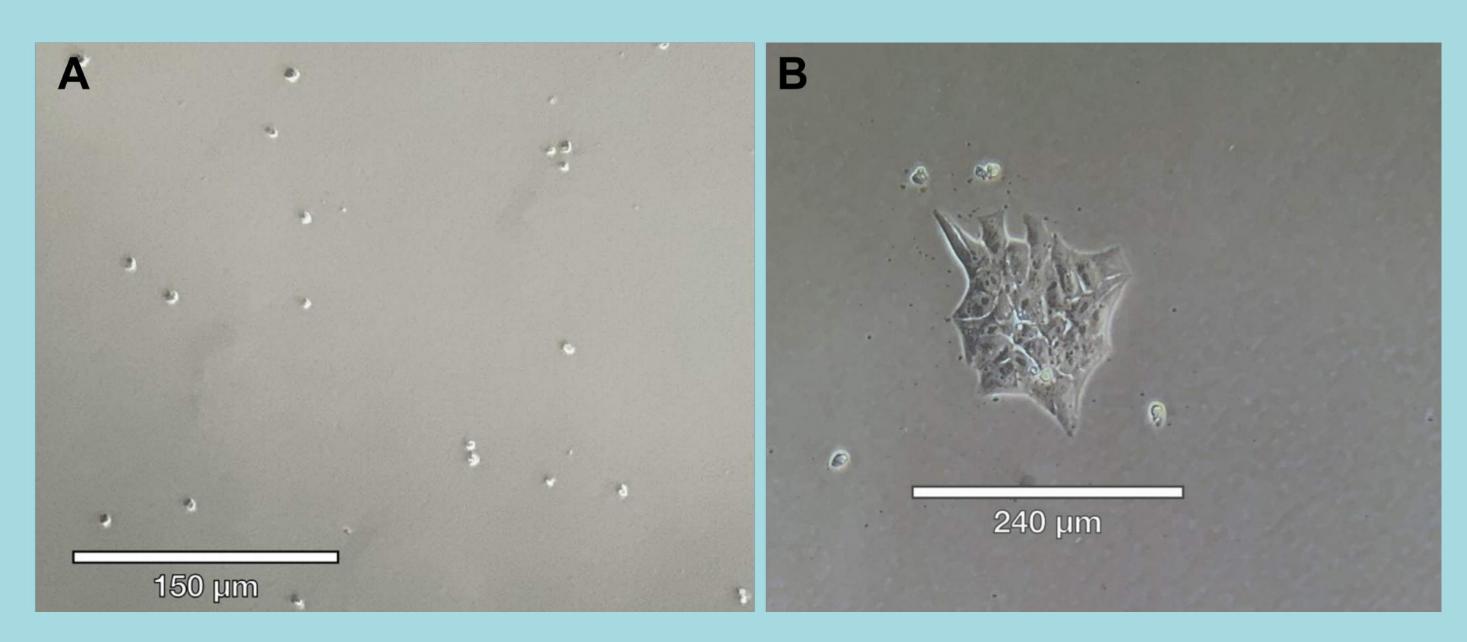


Figure 4. PSCs sorted for clonal expansion: Panel (A) is representative of a WOLF N1 'focusing well' initially plated at higher than one-cell-per-well density to aid microscopy - cell adherence shown 24 hours post-seed. Panel (B) shows archetypal PSC colony emerging (typically 16-32 cells), 5 days after seeding.

## Workflow

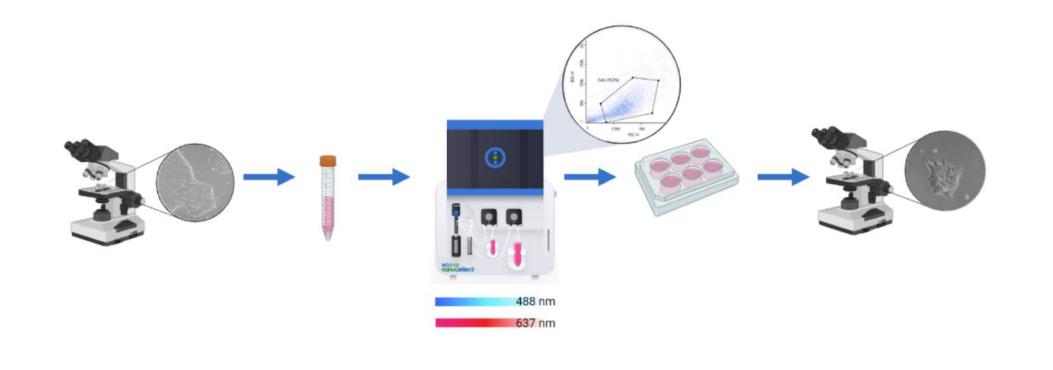


Figure 5. A simplified workflow: The experiment began with general PSC maintenance in HiDef-B8cGMP complete medium and visual confirmation of healthy cells. A single-cell suspension was obtained using defined cGMP-compliant TrypLE Select. PSC suspension was separated into various staining cell tubes and applied to the WOLF G2. Appropriate fluorochromes were chosen to be compatible with the WOLF G2 637 nm and 488 nm LASERs. Double-stained single cell suspension is re-plated into pre-coated plates and imaged on subsequent days for cell health and growth.

# HiDef-B8-cGMP Medium

